


Vitreia® Enterprise Suite PET/CT Viewer

The PET/CT Viewer launches directly from Vitrea® Enterprise Suite, providing enterprise-wide access to PET/CT review tools.

Key features include the following:

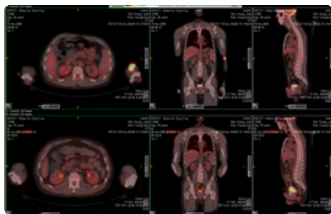

- Specialized PET/CT Layouts
- SUV-based intensity adjustments (W/L presets)
- Volumetric SUV computations
- PET color overlay with blending control
- Ability to load and compare up to 3 time points

PET/CT Viewer automatically selects the paired PET/CT series to load in your viewer.

If your DICOM image data doesn't include key patient information you need to measure SUV, the PET/CT Viewer allows you to easily enter the required data.

Select from four layouts and numerous presets (W/L and color table) to quickly initialize the case presentation. Change the fusion blend between CT and PET data to aid understanding of function and structure.

Calculate and track SUV over time with the adjustable spherical region that makes analysis of potential tumors easy to measure and track with up to three time-points.

Capture and send key images to your PACS or with referring clinicians.

